

For research use only

Catalog NO.: FY-EM3737 Mouse GLP-2(Glucagon Like Peptide 2) ELISA Kit

For the quantitative determination of Mouse GLP-2 concentrations in serum, plasma, and other biological fluids.

This package insert must be read in its entirety before using this product. For research use only. Not for use in diagnostic procedures.

PRINCIPLE OF THE ASSAY

This assay is based on the competitive binding technique. The micro ELISA plate provided in this kit has been pre-coated with Mouse GLP-2. During the reaction, Mouse GLP-2 in samples or Standard competes with a fixed amount of Mouse GLP-2 on the solid phase supporter for sites on the Biotinylated Detection Ab specific to Mouse GLP-2. Following a wash to remove excess conjugate and unbound sample. Then Avidin conjugated to Horseradish Peroxidase (**HRP**) are added to each microplate well and incubated. Substrate solution is added to the wells. The enzyme-substrate reaction is terminated by the addition of stop solution and the color change is measured spectrophotometrically at a wavelength of 450±3nm. The intensity of the color is inversely proportional to the concentration of GLP-2 in the sample.

PERFORMANCE

★ Sensitivity:46.88 pg/mL.

★ Detection Range:78.13-5000pg/mL.

 \star No significant cross-reactivity or interference between Mouse GLP-2 and analogues was observed.

★ Repeatability: Coefficient of variation is 10%.

PRECAUTIONS BEFORE THE EXPERIMENT

1. The detection range of the kit is different from the concentration range of the substance to be measured in the sample. It is recommended to estimate the concentration by looking up the reference before the experiment, and determine the actual concentration by pre-experiment. If the concentration is too high or too low, dilute or concentrate the sample appropriately.

2. The newly opened ELISA plate may contain water fog-like substances, which is a normal phenomenon and will not have any impact on the experimental results.

3. Standard working fluid, biotinized antibody working fluid and HRP Conjugate working fluid should not be reused.

4. The micro plate reader requires a filter with a wavelength of 450 ± 3 nm and a detector capable of detecting 450 ± 3 nm, and the optical density value is between 0 and 3.5.

5. Avoid cross-contamination, change the pipette head between each standard concentration addition, between sample addition, and between reagent addition. In addition, separate storage bottles should be used for each reagent.

6. In order to obtain accurate experimental results, it is necessary to ensure that the sealing plate membrane seals the enzyme label plate during incubation.

7. The substrate solution should remain colorless and free from light until it is added to the plate. When added to the plate, the substrate solution should change from colorless to a gradient of blue.

8. The stop solution should be added to the label plate in the same order as the substrate solution. After the stop solution is added, the color of the solution in the well will change from blue to yellow. Green wells indicate that the stop solution is not sufficiently mixed with the substrate solution. If the color is uneven, shake the plate gently to mix the solution evenly.

REAGENTS&MATERIALS PROVIDED

Item	Specifications	Storage
	96T: 8 wells ×12 strips	
Micro ELISA Plate (Dismountable)	48T: 8 wells ×6 strips	
	24T: 8 wells ×3 strips	
Reference Standard	96T: 2 vials	
	48T: 1 vial	-20°C, 6 months
	24T: 1 vial	
	96T: 1 vial, 60μL	
Concentrated Biotinylated Detection Ab (100×)	48T: 1 vial, 60μL	
A0 (100^)	24T: 1 vial, 60µL	
	96T: 1 vial, 120μL	
Concentrated HRP Conjugate (100×)	48T: 1 vial, 60μL	-20°C, 6 months
	24T: 1 vial, 60µL	
Reference Standard & Sample	1 . 1 . 0	
Diluent	1 vial, 20 mL	
Biotinylated Detection Ab Diluent	1 vial, 13 mL	2-8°C, 6 months
HRP Conjugate Diluent	1 vial, 13mL	
Concentrated Wash Buffer (25×)	1 vial, 30 mL	
Substrate Reagent	1 vial, 10 mL	2-8°C (Protect from light)
Stop Solution	1 vial, 10 mL	2-8°C, 6 months
Plate Sealer	5 pieces	
Product Description	1 copy	RT
Certificate of Analysis	1 copy	

OTHER SUPPLIES REQUIRED

Microplate reader with 450nm wavelength filter Incubator capable of maintaining 37°C Single or multi-channel pipettes with high precision disposable pipette tips EP tubes Container for Wash Solution Squirt bottle, manifold dispenser Deionized or distilled water Absorbent paper for blotting the microplate

SAMPLE COLLECTION&STORAGE

Serum: Allow blood samples to clot for 2 hours at room temperature or overnight at 2-8°C before centrifugation for 15 min at $2000 \times g$ at 2-8°C. Supernatants should be taken for assay testing immediately or stored at -20°C or -80°C for later use. Avoid repeated freeze-thaw cycles.

Plasma: Collect plasma using EDTA or heparin as an anticoagulant. Centrifuge samples for 15 min at 2000×g at 2-8°C within 30 min of collection. Supernatants should be taken for assay testing immediately or stored at -20°C or -80°C for later use. Avoid repeated freeze-thaw cycles.

Cell culture supernatant : Remove particulates by centrifugation and assay immediately or aliquot and store samples at \leq -20 °C. Avoid repeated freeze-thaw cycles.

Other Types of samples : Please check our official website for the collection and processing method(<u>https://feiyuebio.com/index/e31/66.html</u>).

REAGENT PREPARATION

- 1. Allow all reagents to return to room temperature (18-25°C) before use.
- 2. Wash Buffer: Dilute 20 mL of Concentrated Wash Buffer with 480 mL of deionized or distilled water to prepare 500 mL of Wash Buffer.

3. Standard working solution:

①Centrifuge the Reference Standard at 10000×g for 1 min.

②Add 1mL of Reference Standard & Sample Diluent, let it stand for 10min and then mix gently, that is 5000pg/mL of standard working solution.

③Prepare 7 EP tubes, add 500μ L standard & sample diluent to each tube, take 500μ L from 5000pg/mL to make 2500pg/mL, and so on, the last tube is blank (no more liquid needs to be drawn from the penultimate tube). Use it now!



- 4. **Biotinylated Detection Ab working solution**: Calculate the amount of Biotinylated Detection Ab required (50μ L/ well). The actual preparation should be 100-200 μ L. Centrifuge the concentration tube before use and dilute 100× concentrated Detection Ab to 1× working solution (e.g. 10μ L + 990 μ L Biotinylated Detection Ab Diluent).
- 5.**HRP Conjugate working solution**: Calculate the amount of HRP Conjugate required (100μ L/ well). The actual preparation should be $100-200\mu$ L. Centrifuge the concentration tube before use and dilute $100\times$ concentrated HRP Conjugate to $1\times$ working solution (e.g. 10μ L + 990 μ L HRP Conjugate diluent).

ASSAY PROTOCOL

1. Determine wells for diluted standard, blank and sample. Add 50μ L each dilution of standard, blank and sample into the appropriate wells (It is recommended that all samples and standards be assayed in duplicate). Immediately add 50μ L of Biotinylated Detection Ab working solution to each well. Cover the plate with a new sealer. Incubate for 30 min at 37°C. Note: solutions should be added to the bottom of the micro ELISA plate well, avoid touching the inside wall and causing foaming as much as possible.

2. Decant the solution from each well add 300μ L of wash buffer to each well. Soak for 0.5 min and aspirate or decant the solution from each well and pat it dry against clean absorbent paper. Repeat this wash step 3 times. Note: a microplate washer can be used in this step and other wash steps. Make the tested strips in use immediately after the wash step. Do not allow wells to be dry.

3. Add 100µL of HRP Conjugate working solution to each well. Cover the plate with a new sealer. Incubate for 30 min at 37°C.


4. Decant the solution from each well, repeat the wash process for **3 times** as conducted in step 2.

5. Add 100µL of Substrate Reagent to each well. Cover the plate with a new sealer. Incubate for about 15 min at 37°C. Protect the plate from light. Note: the reaction time can be shortened or extended according to the actual color change, but not more than 30 min. Preheat the Microplate Reader for about 15 min before OD measurement.

6. Add 50μ L of Stop Solution to each well. Note: adding the stop solution should be done in the same order as the substrate solution.

7. Determine the optical density (OD value) of each well at once with a micro plate reader set to **450 nm**.

SUMMARY FOR PROCEDURE

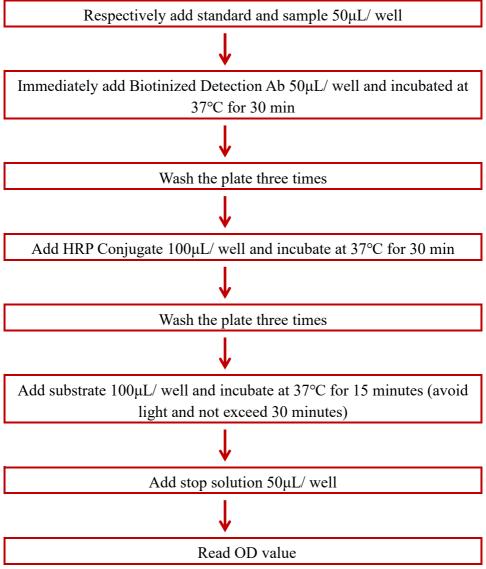

EXAMPLE DATA

For each test, a standard curve must be set for each plate. The standard curves below are for example purposes only.

pg/mL	Corrected OD	Typtical Standard Curve
0	2.552	
78.13	2.075	10
156.25	1.758	and a second
312.5	1.306	
625	0.881	ō _{0.1}
1250	0.549	0.01
2500	0.329	10 100 1000 10000 100000 Mouse GLP-2 Concentration(pg/mL)
5000	0.191	

CALCULATION OF RESULTS

Calculate the average OD value of the standard and sample

With the standard product concentration as the horizontal coordinate and OD value as the vertical coordinate, the standard curve is formed by using the nonlinear four parameter logistic(4-PL) and the sample concentration is calculated by the curve.Please refer to our official website for the specific fitting method(https://feiyuebio.com/index/e31/65.html).

If samples have been diluted, the concentration read from the standard curve must be multiplied by the dilution factor.

SENSITIVITY

The minimum detection concentration of the Mouse GLP-2 detected by this kit is usually less than the 46.88 pg/mL. The sensitivity was determined by adding two standard deviations to the mean OD value of twenty zero standard replicates and calculating the corresponding concentration.

PRECISION

Mean coefficient of variation for Intra-Assay and Inter-Assay: 3 samples with low, middle and high level concentration were tested for repeat multiple times, respectively. The results showed that the coefficient of variation of the kits was less than 10%, which met the precision quality control standard.

SPECIFICITY

This method has high sensitivity and specificity for the detection of Mouse GLP-2. No significant cross-reactivity was observed between the Mouse GLP-2 and its analogues.

RECOVERY

The known concentration of Mouse GLP-2 was added to different samples respectively, and the recovery experiment was conducted. The results showed that the recovery range and average recovery rate of the kit were 80-120%, which met the recovery quality control standard.

LINEARITY

Dilute linear experiments were performed on samples with Mouse GLP-2 to evaluate the linearity of the kit. The results showed that the kit linear range (%) was 80-120%, which met the linear quality control standard.

DECLARATION

1. Limited by current conditions and scientific technology, we can't conduct comprehensive identification and analysis of all the raw material provided. There might be some qualitative and technical risks for users using the kit.

2. This assay is designed to eliminate interference by factors present in biological samples. Until all factors have been tested in the ELISA immunoassay, the possibility of interference cannot be excluded.

3. The final experimental results will be closely related to the validity of products, operational skills of the operators, the experimental environments and so on We are only responsible for the kit itself, but not for the samples consumed during the assay. The users should calculate the possible amount of the samples used in the whole test. Please reserve sufficient samples in advance.

4. To get the best results, please only use the reagents supplied by the manufacturer and strictly comply with the instructions.

5. Incorrect results may occur because of incorrect operations during the reagents preparation and loading, as well as incorrect parameter settings of the Micro plate reader. Please read the instructions carefully and adjust the instrument prior to the experiment.

6. Even the same operator might get different results in two separate experiments. In order to get reproducible results, the operation of every step in the assay should be controlled.

7. Every kit has strictly passed QC test. However, results from end users might be inconsistent with our data due to some variables such as transportation conditions, different lab equipment, and so on. Intra assay variance among kits from different batches might arise from the above reasons too.

8. Kits from different manufacturers or other methods for testing the same analyte could bring out inconsistent results, since we haven't compared our products with those from other manufacturers.

9. The kit is designed for research use only, we will not be responsible for any issues if the kit is applied in clinical diagnosis or any other related procedures.

PLATE LAYOUT

Use this plate layout to record standards and samples assayed.



If you have any questions, please contact us. Email: <u>info@feiyuebio.com</u> <u>tech@feiyuebio.com</u>